Grid Pad 15

Created by

<u>Smart</u>box

Introducing Grid Pad 15

A fully featured communication aid designed for independence. It works with all access methods and features a 15.6" display and rear-facing second screen.







- **Communication**
- **T** Environment control
- **Computer control**
- Accessible apps

Purpose built for independence

With a large 15.6" display and full colour second screen, Grid Pad 15 is a flexible solution for people with complex access and communication needs. A long 14-hour battery life, rugged design and multiple access options ensure you always have your voice.















Touch

Pointer

Switch

Eye gaze

- ✓ High spec computer for fast communication
- √ 15.6" HD tough antiglare display
- Antimicrobial casing and screen
- √ 4.6" rear-facing second screen
- Built for all access methods
- ✓ Over 14 hours battery life (9 hours with eye gaze)
- High quality amplified sound
- Remote power button
- ▼ Built-in radio and infrared Environment Control
- ▼ Tactile, water-resistant buttons
- Flexible stand and mounting options

Key features



Tough touchscreen display

Grid Pad 15 has a 15.6" HD touchscreen, with a hardened glass display. The antiglare finish helps to reduce reflections and makes it easier to see the screen when you are out and about.

Powerful device

Grid Pad 15 has a fast processor, perfect for communicating quickly and running programs using Computer Control.

All day battery

Grid Pad 15 is designed to help you communicate all day long. The large battery provides up to 14 hours continuous use (9 hours with eye gaze).



Antimicrobial casing and screen

Antimicrobial Sanitized® technology within the case and screen inhibits the growth of bacteria and viruses.

Switch and USB devices

There are two built-in ports for switches, and two USB ports to use with any other type of alternative access, such as a joystick, trackball or HeadMouse Nano.

Eye gaze ready

Grid Pad 15

Turn your Grid Pad into an integrated eye gaze device with our purpose-built camera unit. Available with Lumin-i or Alea eye gaze technology.



We also have bespoke mounting plates for use with other eye trackers.



Full colour second screen

The 4.6" second screen shows your message on the rear of the device, improving communication in a range of situations. Animated dots let others know not to interrupt while you compose your message.

Share your messages privately rather than projecting audio to everyone around and control when others can see it with one click.





Antimicrobial casing

Grid Pad 15 is built with an antimicrobial casing and screen to inhibit the growth of bacteria and the spread of viruses, helping to keep your device cleaner for longer.

The technology has not yet been tested on COVID-19, however it has been successfully tested against similar viruses including Influenza, Avian flu and SARS.

- ✓ Inhibits bacteria and virus growth
- Prevents unpleasant odours
- Cannot be removed by cleaning
- ✓ Sanitized® technology





Accessories





Remote power button

Grid Pad comes with our remote power button, to give you more independence. You can plug any switch into the button and use it to turn the Grid Pad 15 on.




Need a keyguard?

Our latest keyguards use non-slip suction technology to fix in place. Super easy to assemble and clean, no tools required. Available as an optional extra.

Manual, dual mounting plate, screwdriver, power lead, Grid 3 getting started guide, cleaning cloth, optional port labels

Other key features



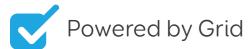

Everything you need for communicating, controlling your environment, message banking, browsing the internet, using social media and much more!



The integrated speakers are engineered to produce high quality sound and volume, to help get your message across effectively wherever you are.



Easily mount your device to a floor stand with the Grid Pad dual mounting plate, compatible with Daessy and Rehadapt mounting systems.



Control everything from your lights to your TV.
The infrared learning window is located on the front for quick and easy programming.





Operating system	Windows 10 Pro
Software	Grid 3
Processor	Intel® Core™ i5-7Y54 Processor (1.2GHz, 3.2GHz Burst)
Memory / RAM	8GB
Hard drive	SSD SATA III, 256GB
Display	15.6" 1920 x 1080 with toughened glass and anti-glare finish
Second screen	4.6" 800 x 300 with toughened glass and anti-glare finish
Weight	3kg 6.6 lbs
Battery life	14 hours typical use 9 hours with eye gaze
Charge time	10 hours

Tablet camera	2MP front, 5MP rear
Dimensions	408 x 259 x 42 mm 21 x 10 x 1.7"
Eye gaze	Optional Grid Pad camera unit (Alea or Lumin-i) Eye gaze mounting plates also available
Access	2 x 3.5mm switch ports, 3.5mm headphone / mic combo port, touch screen, magnetic Keyguard compatible, 1 x USB 2.0 port, 1 x USB 3.0 port (Eye gaze uses the USB 3.0 port)
Environment control	GEWA Infrared, EasyWave / Z-Wave radio
Connectivity	Dual band wireless (AC 8265), Bluetooth 4.2 LE
Mounting plate	Daessy and Rehadapt
Warranty	3 years SmartCare (UK only)
Other features	Antimicrobial casing and screen, built in stand and remote power button

thinksmartbox.com